

ß Gal-Juice/ ß Gal-Juice PLUS (1 033 11/ 1 033 12)

Components included:

β-Gal-Juice/plus 50ml

Chemiluminescent Substrate for measurement of ß-Galactosidase

Enzyme.

Store at +4°C.

Triggering Reagent 50ml

Store at +4°C

Preparation of Cell Lysates:

For cell lysis, we suggest the follow buffers:

• **Animal cells**: 100 mM potassium phosphate (pH 7,8), 0,2% Triton X100

• Yeast cells: 0,5 M sodium phosphate (pH 7,1) 50 mM KCl, 5 mM MgSO₄

Standard Protocol for Detection of β-Galactosidase in Microtiter Plate Luminometer

Equilibrate at room temperature for 30 minutes before use.

Best results for chemiluminescent detection of β -Galactosidase enzyme or reporter assays can be obtained from 30 minutes to 60 minutes incubation of β -Gal Juice with β -Galactosidase enzyme and read the plate or tube within 2 to 10 minutes after triggering the reaction mixture by accelerator or enhancer.

- Add 5 to 10 μl of diluted β-Galactosidase enzyme or cell extract to microplate wells.
- Add 50 to 100 μ l of β -Gal Juice and incubate for 30 to 60 minutes at room temperature.
- Add 50 to 100 μl of triggering reagent and shake it for 10 seconds.
- Place the microtiter plate in luminometer and start reading.
- Best results are obtained between 2 to 10 minutes.



ß Gal-Juice/ ß Gal-Juice PLUS (1 033 11/ 1 033 12)

Components included:

β-Gal-Juice/plus 50ml

Chemiluminescent Substrate for measurement of ß-Galactosidase

Enzyme.

Store at +4°C.

Triggering Reagent 50ml

Store at +4°C

Preparation of Cell Lysates:

For cell lysis, we suggest the follow buffers:

• **Animal cells**: 100 mM potassium phosphate (pH 7,8), 0,2% Triton X100

• Yeast cells: 0,5 M sodium phosphate (pH 7,1) 50 mM KCl, 5 mM MgSO₄

Standard Protocol for Detection of β-Galactosidase in Microtiter Plate Luminometer


Equilibrate at room temperature for 30 minutes before use.

Best results for chemiluminescent detection of β -Galactosidase enzyme or reporter assays can be obtained from 30 minutes to 60 minutes incubation of β -Gal Juice with β -Galactosidase enzyme and read the plate or tube within 2 to 10 minutes after triggering the reaction mixture by accelerator or enhancer.

- Add 5 to 10 μl of diluted β-Galactosidase enzyme or cell extract to microplate wells.
- Add 50 to 100 μ l of β -Gal Juice and incubate for 30 to 60 minutes at room temperature.
- Add 50 to 100 μl of triggering reagent and shake it for 10 seconds.
- Place the microtiter plate in luminometer and start reading.
- Best results are obtained between 2 to 10 minutes.




ß Gal-Juice/ß Gal-Juice PLUS (1 033 11/ 1 033 12)

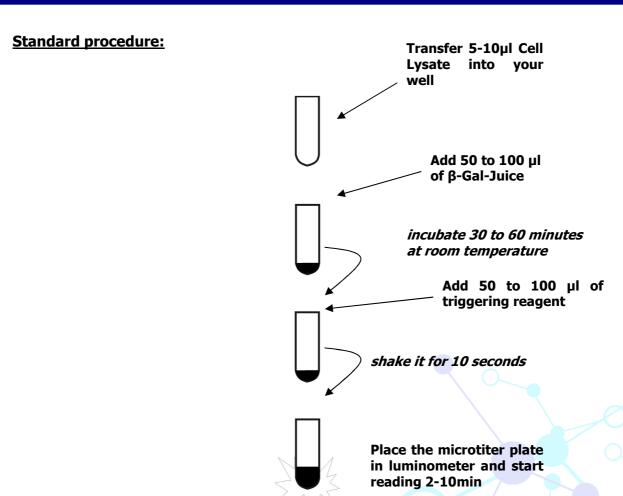

Optional: parallel measurement in Duo- or Triple-Well System (Firefly/ Renilla/ Gaussia/ AP/ βGal)

Split the cell lysate and transfer to different tubes for measuring several enzymes in one sample.

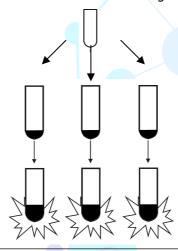


PJK GmbH

Elsaesserstr. 51 - 66271 Kleinblittersdorf - Germany

Phone: +49 (0)6805-91 20 22 - Fax: +49 (0)6805- 91 20 24




ß Gal-Juice/ß Gal-Juice PLUS (1 033 11/ 1 033 12)

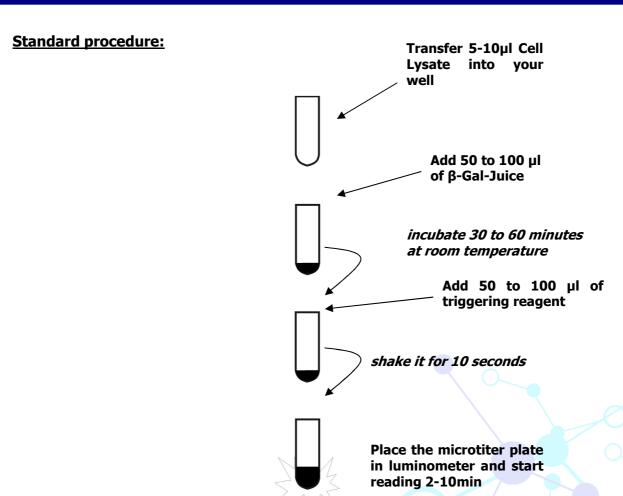

Optional: parallel measurement in Duo- or Triple-Well System (Firefly/ Renilla/ Gaussia/ AP/ βGal)

Split the cell lysate and transfer to different tubes for measuring several enzymes in one sample.

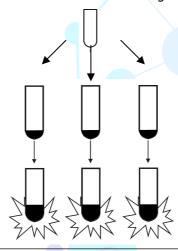

PJK GmbH

Elsaesserstr. 51 - 66271 Kleinblittersdorf - Germany

Phone: +49 (0)6805-91 20 22 - Fax: +49 (0)6805- 91 20 24